

*Brighter
Your Smile ...!*



KORdENT

Precision makes innovation



■ Restorative Materials

- K. Composite Resin
- K. Composite Resin Kit
 - Kit 1, Kit 2
- K. Flow Resin
- K.Temp Flow
- K.Tem
- K. Temp Crown
- K. Temp Crown Kit
- K.flow for Implant
- K.Lay

■ Etchant & Adhesives

- K. Etch-37
- K. Bond
- K. Temp Cement
- K. Temp Cement NE
- K. Temp Cement EZ
- K. Temp Cement Implant

■ Endodontic Materials

- K. Sealer
- K. Bio ceramic sealer
- K. Pex
- K. Paste
- K. MTA

■ Impression Materials

- K. Light Body
- K. Putty Set
- K. Alginate
- K. Tray Cleaner
- K. Sil



■ Core Materials, Liner EDTA

- K. Core
- K. Base Liner
- K. EDTA

■ Prevention & protection Pumice

- K. Fluoride Varnish
- K. Varnish
- K. Fluoride Gel
- Flavors : strawberry, orange, peach
- K. Fluoride Gel Tray
Small, Medium, Large
- K. Pumice
- Pumice Paste (with fluoride)
- Pumice Paste
(without fluoride)

■ Tooth Whitening Technique

- K. Bleaching
- Kit 1, Kit 2

■ Another Products

- K. Cord (Retraction Cord)
- K. Oil
(lubricant Handpieces)
- K. Dam



Precision makes innovation

Kordent is a provider of full range of dental products for global partners with the proven track of records.

Kordent was founded on the principle to provide utmost satisfaction for customers by offering the best products at the best value. It is also **Kordent's** deepest intention to build up the healthy

dental care for people around the world.

The increasingly tough competition in the world market has consistently stimulated **Kordent** to raise its standards to the highest quality to maintain the best value for price, and to stay at the forefront of the latest up-to-date technology to satisfy the high-end product expectancy.

Kordent's drive to upgrade and maintain the high quality is a continuous process. And it wishes to share its drive with the global partners who are set to reap the benefits in their markets with what **Kordent** offers through its Denu brand.

Kordent has always upheld the value to satisfy its customers as the number one priority. And that value will continue to be the core of **Kordent** principle.

K. Composite Resin / Composite Resin Kit

- K. Composite Resin

Light curing universal composite resin with nano hybrid filler

K. Composite Resin is water insoluble with the optimal biocompatibility and provides the high polymerization rate, low shrinkage with enhanced flexural rigidity. Stickiness and slumping is minimized to allow excellent usability for the application of all anterior and posterior teeth restorations (Class I~V)

Shade: A0, A1, A2, A3, A3.5, B1, B2, B3, U0

Features

- High strength and excellent water insoluble with optimum nano filler
 - Compressive strength 396 Mpa, Flexural strength 175 Mpa, Hardness 95 Mpa
- Excellent applicability (No stickiness and easy application)
- Manifested clearly on X-ray scan, Highly polishable
- Applicable for all anterior and posterior teeth restoratic(Class I~V)
- Excellent color match with natural tooth

Suggested Applications

- Restoration for all anterior and posterior teeth sections(Class)-V
- Restoration to fractured anterior teeth
- Replacement of amalgam for esthetic purpose
- Core build-up

Package

- 1 x 4g syringe/pack

Clinical Procedures



Preparation



Etching



Clean and dry




Bonding



Lightening



Filling

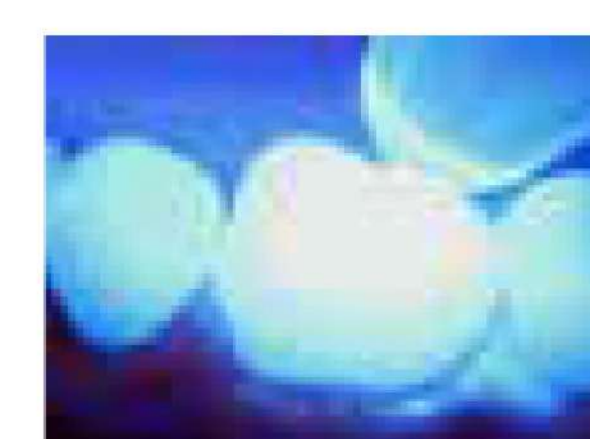

Lightening



Polishing

- K. Composite Resin Kit

Kit 1

Composition

- K. Composite Resin - 8 syringes
- K. Etch-37-1 syringe, tips
- K. Bond-2 bottles, 15 micro applicators
- * Selectable composite resin shades
- (Basic shades: A0, A1, A2, A3, A3.5, B1, B2, B3, U0)

Kit 2

Composition

- K. Composite Resin - 2 syringes
- K. Etch-37 - 1 syringe, 2 tips
- K. Bond- 1 bottle



Precision makes innovation

Restorative Materials

High Quality With Kordent

K. Flow Resin / K. Temp Flow

- K. Flow Resin

Light cured flowable composite resin with nano hybrid filler

The flowability ensures easy applications that require details and also easy to use on tooth bridge or tooth barriers. It is also used as Base and Liner for restorations as well as for tooth bridges or tooth barriers

Shade: A0, A1, A2, A3, A3.5, B1, B2, B3, U0

Features

- Excellent water insoluble material with optimal nano hybrid filler
 - Compressive strength 340 Mpa, Flexural strength 150 Mpa, Hardness 65 Mpa
- Excellent flow into the dental cavities by superb flowability
- Low polymerization shrinkage, Available in adhering orthodontic bracket
- Aesthetic is well considered by matching the near-natural tooth

Suggested Applications

- Cavity liner for direct restoration, Small and deep minute gap
- Sealing of pit and fissure, Base and Liner of restoratives
- Block-out of undercut

Package

- 2 x 2g syringes, disposable tips / pack

- K. Temp Flow

Light curing semi-gel type temporary composite resin

Used as temporary filling material to isolate prepared dental cavity from the extreme oral conditions before dental prosthetic installation. Semi-gel type to be directly discharged to the dental cavity through disposable tip for easy application

Shade: Blue, Light Blue, Yellow

Features

- No easily deviation till dental prosthetic installation
- Semi-gel type temporary restoration material
- Easy to remove, No debris left after removal
- Excellent sealing ability, Adhesive to the side of dental cavities
- Easy to check debris after removal with visible color (blue, yellow)

Suggested Applications

- Temporary restoration for inlay, onlay cavity
- Temporary filling and sealing, Temporary filling for implant abutment
- Damp-blocking temporary sealing,

Package

- 5 x 2g syringes, disposable tips / pack



Precision makes innovation

Restorative Materials

High Quality With Kordent

K. Tem / K .Temp Crown

- K. Tem

Light cured resin for making a temporary tooth for permanent prosthesis or implant

Easy & fast fabrication of temporary tooth by hand,
reducing the material loss

Shade : A2, A3

Features

- Light cured resin for fabrication of temporary tooth
- Fast fabrication (within 3 minutes) of temporary tooth by hand, reducing the material loss
- Easy to use, Syringe type

Package

- 5 x 12ml (14g) syringes / pack
[One syringe can make approximately 20 temporary teeth or 30 temporary anterior teeth]

- K. Temp Crown

Light cured paste-type resin for making a temporary crown

Shade : A1, A3

Features

- Enables to save the materials and time for fabrication of temporary crown

Package

- 2 x 2g syringes, tips

- K. Crown Kit

Temp Crown Kit for fabrication of temporary crown

Composition

- K. Temp Crown : 10 x 2g syringes, tips
- K.Trans Sil : 2 x 50ml (25ml base, 25ml catalyst) cartridges
- K. Vaseline : 1 x 5ml syringe, tip

Apply Vaseline on the tooth and set Crown Resin after making a shape with Transparent Silicone and light cure Transparent Silicone from the outside



Precision makes innovation

Restorative Materials

High Quality With Kordent

K. flow for Implant / K. Lay

- K.flow for Implant

Flowable filling material for implant abutment

Features

- Perfect sealing to prevent screw loosening
- Excellent homogenous for perfect sealing
- No shrinkage, No expansion
- Flexible enough to remove easily
- Easy to apply, No cracking, No foul odor

Suggested Applications

- Temporary filling material for **implant abutment**

Package

- 5 x 1.2ml syringes, disposable tips

- K.Lay

Light cured semi-flowable block-out material

Features

- Low viscosity for easy flow
- Fast working , Easy to remove
- Blue color for easy visibility

Suggested Applications

- Blocks out undercuts on dies and fills in voids
- Provides reservoir space for bleaching individual trays
- Blocks out the part to be protected before sanding

Package

- 5 x 1.2ml syringes, disposable tips

K. Etch-37 / K. Bond

- K. Etch-37

Semi-gel type etching agent with phosphoric acid

Features

- Optimal viscosity for precise placement and excellent control
- Excellent etching of dentin and enamel
- Provides excellent etching depth
- Perfect removal and no debris left after etching
- No leak from the syringe after use

Suggested Applications

- Etching for bonding composites, sealants or adhesives

Package

- 5 x 3ml syringes, disposable tips

Clinical Procedures



Preparation



Etching



Clean and dry

- K. Bond

One-step for priming and bonding 5th generation bond

Features



- Light cured single component dentin and enamel bonding agent
- 5th generation bond (Use Denu Bond after etching and then the resin)
- High bond strength, Long-lasting bond
- Available in bottle delivery

Suggested Applications

- All direct placements for composite resin restoration
- Bonding of composite to dentin, enamel and cementum

Package

- 2 x 5ml(5g) bottles

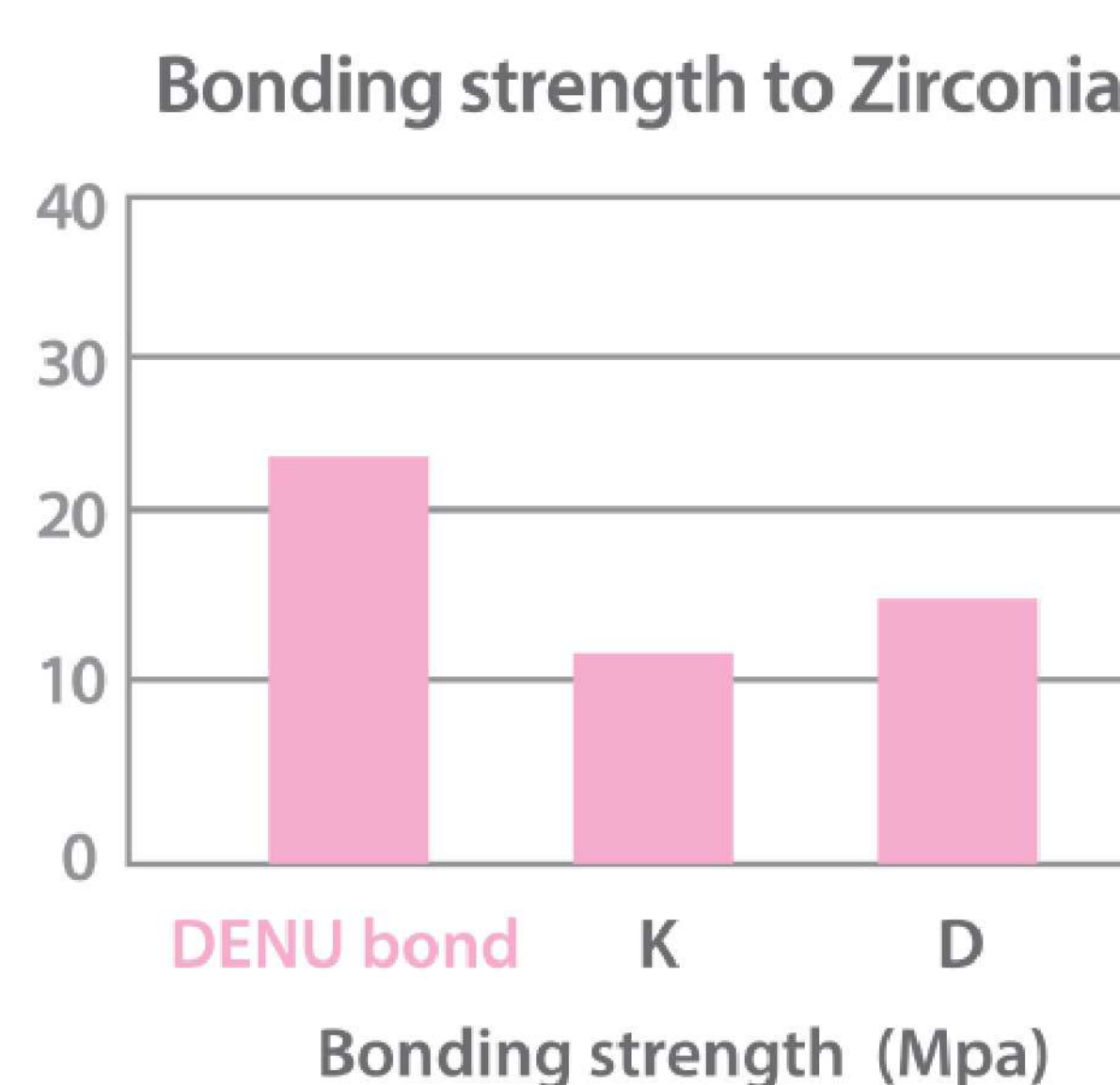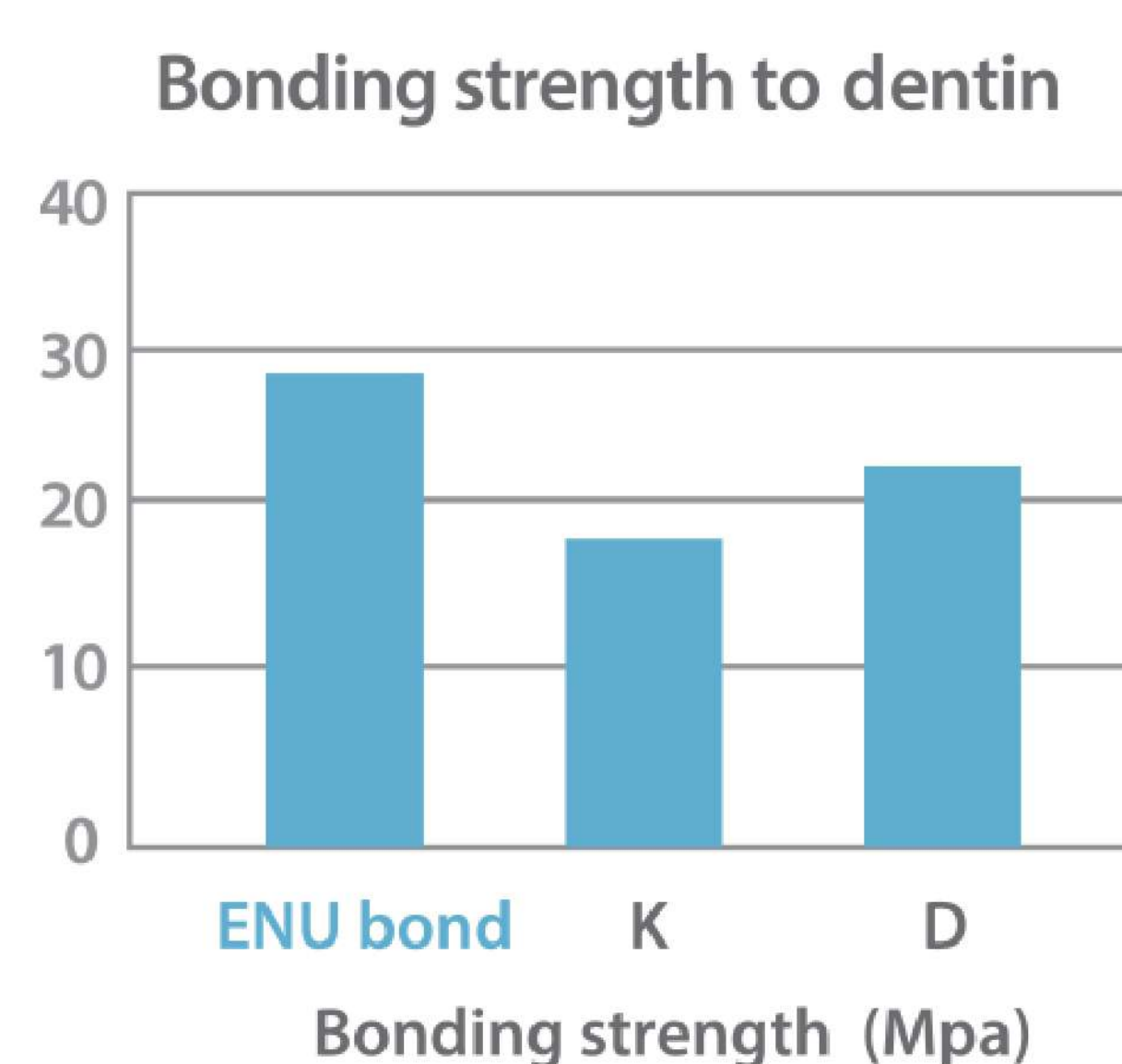

K. Temp Cement / K. Temp Cement NE

- K. Temp Cement

Temporary cement solvent, Removable up to 0.3mm thick


Features

- Perfect removal of temporary cement
- Applicable for eugenol or non-eugenol temporary cement
- No damage to prosthesis
- No contamination by debris

Suggested Applications

- Removal of temporary cement

Leave the prosthesis such as crowns to soak in Denu Temp Clean for 10 minutes and clean it with running water for 10 seconds.

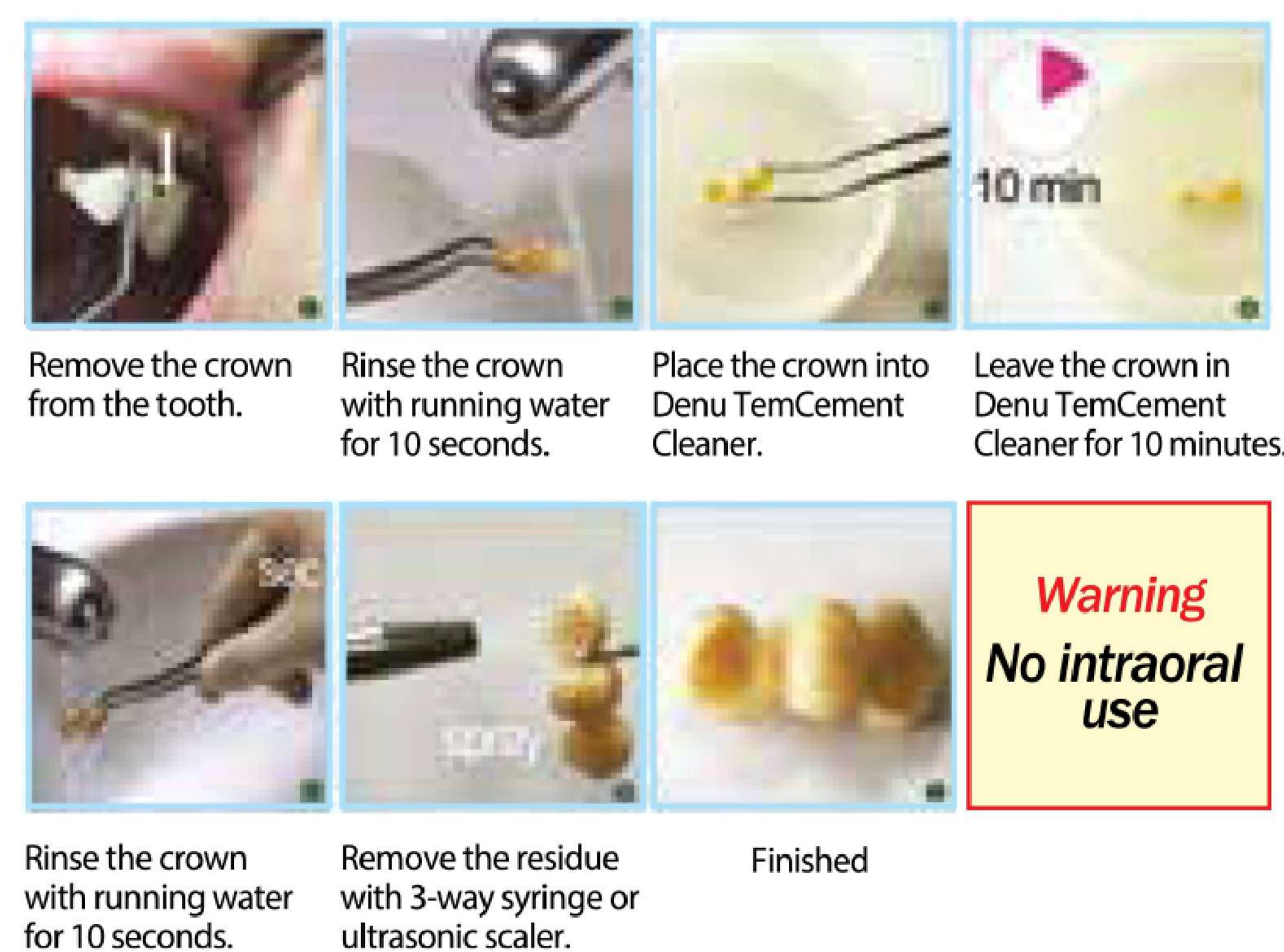

Package

- 1 x 500ml bottle
- 1 x 200ml bottle

- K. Temp Cement NE

Non-eugenol temporary cement

Features

- Good adhesive bonding strength
- Compatible with acrylic resin and plastic provisional materials
- With attachment of eco tip, minimum amount of Denu Temp Cement is used
- Due to fast curing in oral cavity, operation time is reduced
- No irritation of soft tissue

Suggested Applications

- Temporary cementing of temporary ceramic, metal or composite inlay, onlay, crown, bridge, post and screw

Easy Removal

Package

- 2 x 5ml(13g) syringes, eco tips
- Mixing time: 10-15 seconds, Operation time : within 30 seconds
- Setting time in mouth: 3.30 minutes



Precision makes innovation

Etchant & Adhesives

High Quality With Kordent

K. Temp Cement EZ / K. Temp Cement Implant

- K. Temp Cement EZ

Non-eugenol temporary cement

Features

- Easy to remove residual cement on inner surface of prosthetic appliance or tooth surface
- Mixing time: 10-15 seconds, Operation time: within 40 seconds
Setting time in mouth: 3 minutes

Package

- 2 x 5ml(13g) syringes, eco tips

Dual cure type Radiopaque

- K. Temp Cement Implant

Non-eugenol temporary resin cement for implant

Features

- Good adhesive bonding strength
- Dual-cured (light-cured, self-cured)
- Compatible with acrylic resin and plastic provisional materials

Suggested Applications

- Temporary cementing of temporary prosthetic appliances for implant

Excellent cleaning within 10 minutes

Package

- 1 x 5ml(13g) syringes, eco tips
- Operation time: within 1.30 minutes
Setting time in mouth: 4 minutes



Precision makes innovation

Endodontic Materials

High Quality With Kordent

K. Sealer / K. Bio ceramic sealer

- K. Sealer

Root Canal Sealing Material - Epoxy Resin Based

Features

- The ideal root sealer with Low Shrinkage reduced gapping between the sealer and canal wall.
- Low solubility makes it resistant to breaking down over time.
- 1:1 paste to paste manual mixing system.
- Fast and Easy Preparation and Less Waste.
- Excellent flow properties allow it to reach and seal lateral canals,
- Low risk of cytotoxicity with short setting time.
- High radiopacity
- Lightly coated Gutta Percha Points with a good film thickness.

- K. Bio ceramic sealer

k. BioCeramic sealer Sealing Material

Features

DIA-ROOTT™ BIO Sealer is a bioceramic MTA-based (Calcium aluminate-based), non-shrinking root canal sealer. It is used for the permanent obturation of root canals and suitable for all obturation techniques involving Gutta-Percha. is a bioceramic MTA-based (Calcium aluminate-based), non-shrinking root canal sealer. It is used for the permanent obturation of root canals and suitable for all obturation techniques involving Gutta-Percha.

Package

- 1 syringe x 2g
- 20 Disposable tips
- 1 Silicon cap



Precision makes innovation

Endodontic Materials

High Quality With Kordent

K. Pex / K. Paste / K. Mat

- K. Pex

Calcium Hydroxide Paste with Iodoform

Benefits

- Quiets "Hot Abscess"
- Disinfects canal
- Promotes Apexification
- Treats Traumatic Injuries
- Radiopaque
- Quick and simple application
- Antibacterial and Bacterostatic

Indications

- Intracanal medicament
- Apexification -
- Exudation Control
- Periapical Lesions
- Root Resorption
- Temporary root filling
- Perforations
- Under developed Pulpless Teeth

- K. Paste

Calcium Hydroxide Paste with Barium Sulfate

Benefits

- Water Soluble:
- Easy to clean and remove, reduce operative time
- Easy and quick access
- Radiopaque
- Antibacterial

Characteristics

- Premixed Calcium Hydroxide + Barium Sulfate paste for root canal treatment.
- Temporary root filling.
- Root canal filling for primary teeth.
- Use after pulp extirpation treatment.
- Filling on the pulpotomy.
- Temporary Pulp Capping.

- K. Paste

Mineral Trioxide Aggregate

Benefits

- Good Biocompatibility
- Resists washout
- Easy Handling Properties
- Outstanding Sealing Capabilities
- White powder is ideal for Anterior teeth
- Antibacterial Effect (pH 12)

Package

- ITEM #1003-601
- MTA Powder 0.5g
- 1 Mixing Pad
- 1 Spatula



Precision makes innovation

Core Materials, liner EDTA

High Quality With Kordent

K. Core / K. Base Liner / K. EDTA

- K. Core

K. Core flowable core build-up resin

Features

- Dual-cured (light cured, self cured)
- Fast and deep dual-cured (Polymerization depth for 20 seconds : 5.1mm)
- High compressive strength (310 MPa)
- High bond strength to dentin, enamel and post

Package

- 2 x 5ml(10g) syringes, mixing tips, eco tips

Suggested Applications

- Applicable for all types of core build-ups
- Permanent restoration of tooth damage for therapeutic purpose

- K. Base Liner

Light cured cavity liner and pulp-capping material prior to placement of restorative materials

Features

- Pulp-capping material
- Protect a dental pulp for fillings, High strength
- Contains calcium hydroxide in a urethane dimethacrylate base
- Low shrinkage when light-cured

Suggested Applications

- Cavity liner and base material prior to placement of restorative materials
- Pulp-capping

Package

- 2 x 1.2ml syringes, tips

- K. EDTA

Root canal preparation EDTA

Features

- Root canal cleaner with unique formulation of glycol and urea peroxide
- Excellent chelation removes inorganic compound
- Lubrication provides easy instrument handling when root canal is shaped
- Tooth-colored pigment produces no pigmentation
- Perfect removal is possible after application
- Low viscosity of the creamy type provides high flowability and easy spraying into root canal

Suggested Applications

- Cleaning and preparation of root canal
- For expansion of root canal and removal of smear layer

Package

- 2 x 6g syringes, disposable tips



Precision makes innovation

Prevention & protection Pumice

High Quality With Kordent

K. Fluoride Varnish / K. Varnish

- K. Fluoride Varnish

Contains a flavored, xylitol-sweetened, 5% sodium fluoride with TCP and vitamin E in a natural pine resin carrier

Features

- Contains high-enriched fluoride, Long-lasting fluoride release
- Contains tricalcium phosphate(TCP) for continuous and effective fluoride release
- High adsorptive power for fast tooth surface coating
- Long-lasting sealing of the exposed dentinal tubules for relief of hypersensitivity
- Transparency provides natural tooth color for esthetic effect
- Fresh sweet apple flavor for easy access to children

Package

- 10 unit-dose package:
10 unit doses (Varnish and applicator brush are contained in each unit dose)
- 100 unit-dose package:
100 unit doses (Varnish and applicator brush are contained in each unit dose)

Suggested Applications

- Caries prevention and hypersensitive tooth treatment

- K. Varnish

Contains flavored, xylitol-sweetened, 5% sodium fluoride with Vitamin E in a natural resin carrier for effective caries prevention and hypersensitive tooth treatment

Features

Excellent flowability for easy application, Long-lasting fluoride release

- High adsorptive power for fast tooth surface coating, Contains xylitol for effective caries prevention
- Long-lasting sealing of the exposed dentinal tubules for relief of hypersensitivity
- No feeling of irritation in the mouth after application due to the particles on the teeth
- Light yellow color of Denu Varnish indicates the application, but disappears and becomes tooth color in 5 minutes after the application
- Fresh natural orange flavor is added to reduce unpleasant smell and bitter taste of pine resin

Suggested Applications

- Caries prevention and hypersensitive tooth treatment

Package

- 10 unit-dose package:
10 unit doses (Varnish and applicator brush are contained in each unit dose)
- 100 unit-dose package:
100 unit doses (Varnish and applicator brush are contained in each unit dose)



Precision makes innovation

Prevention & protection Pumice

High Quality With Kordent

K. Fluoride Gel / K. Fluoride Gel Tray

- K. Fluoride Gel

Contains fluoride compounds and is used for dental care such as caries prevention by slightly biting a fluoride gel tray with fluoride gel filled

Features

- Natural dental care including caries prevention
- Contains fluoride ion (1.23%), xylitol and vitamin E
- Does not flow down when the gel is filled in the tray
- Effectiveness varies depending on the biting time of the fluoride gel tray filled with the fluoride gel (1 minute: 80% effective, 4 minutes: 100% effective)

Package

- 1 x 500ml bottle

Suggested Applications

- Caries prevention using a fluoride gel tray

- K. Fluoride Gel Tray

Disposable dual-arch tray designed for application of fluoride gel

Features

- Dual arch-shaped tray for lower and upper teeth
- No easy outflow from the tray when the fluoride gel is filled in
- Made of harmless soft rubber to minimize feeling of resistance or irritation
- Curved grip for good grip feeling and tray folding

Suggested Applications

- Application of fluoride gel for caries prevention

Package

- 25 pcs / small box
- 100 pcs / big box (4 small boxes contained)



Precision makes innovation

Prevention & protection Pumice

High Quality With Kordent

K. Pumice

- K. Pumice

Material for grinding and polishing of teeth or prosthesis

Features

- Enables to polish and smooth rough surface of tooth or prosthesis
- Phosphated fluoride (1.23% fluoride ion) is contained
- Optimal viscosity leaves less lumps and provides fast grinding
- Nano-sized natural particles provide fine abrasion
- No pigmentation after application
- Sweet strawberry and natural mint flavors

Package

- 1 x 340g bottle

Suggested Applications

- Grinding and polishing of teeth or prosthesis



Precision makes innovation

Tooth Whitening Technique

High Quality With Kordent

K.Bleaching

- K.Bleaching

Contains water-soluble monomer for gel type bleaching
A single application can generate excellent result

Features

Appropriate viscosity prevents the rundown during treatment
Superb penetration power generates excellent whitening result
Convenient packaging for single treatment
Excellent color match with natural tooth

Package

Whitening Kit 1

- 10 x 1.5ml whitening agent jars, 10 x 3.7g H2O2 bottles
- disposable spatulas, disposable brushes

Package

Whitening Kit 2

- 1 x 1.5ml whitening agent jar, 1 x 3.7g H2O2 bottle
- 2 disposable spatulas, 2 disposable brushes,
2 syringe Dam, case



Precision makes innovation

Another Products

High Quality With Kordent

K. Cord / K. Oil

- K. Cord

is made of 100% cotton, knitted into thousands of tiny loops to form long, interlocking chains. This unique knitted design exerts a gentle, continuous outward force following placement, as the knitted loops seek to open. Optimal tissue displacement occurs in 3–8 minutes. Ultrapak E knitted cord is impregnated with racemic epinephrine hydrochloride.

Features

The original—and smallest—knitted cord

Provides rapid tissue displacement, detailed margins, and quality impressions

Facilitates easy packing and stays placed better than any twisted or braided cord

Compresses upon packing then expands for optimal retraction

Does not entangle in diamond bur

Bright colors facilitate easy location and removal

Available with or without epinephrine*

High carbon steel blade ensures a clean cut

Thin plastic gate prevents cord from falling into the bottle upon cutting

Ruler printed on label provides easy cord measurement

Suggested Applications

Do not use K. Cord or other epinephrine preparations with ferric sulfate solutions, including hemostatics, as a blue/black precipitate will occur.

K. cord with Epinephrine is not available in all countries

- K.Oil

Cleans well with no clotted residue

Features

The product provides excellent cleaning high & low speed handpieces.

Its oil does not clog and is less sticky and does not leave any residue.

Its alcoholic ingredient provides superior disinfection and high volatility making the oil smell disappear in 2-3 seconds after spraying.

The product can be used for all kinds of handpieces including Kavo, W & H, NSK and Bien Air.

Its high compression air in spray container enables highly accurate spraying, which creates minute and uniform coating film on the surface of bearings inside handpiece.

Suggested Applications

1. Shake well before using it.
2. After decomposing air turbine handpiece from tubing coupler, tighten spray nozzle to handpiece.
3. Lubricate for 2 ~ 3 seconds till oil come out.

Package

Spray 520ml×1



Precision makes innovation

Another Products

High Quality With Kordent

K. Dam

- K. Dam

Dam

Light cured gum-protecting material for tooth whitening

Gingival protection resin material out of Light Curing resins to protect gum damage caused by bleaching agents during the tooth whitening process.

Features

Optimal viscosity for easy handling

Short light-curing time (15 seconds)

Adequate flow, easy usability.

Easy to remove after the treatment

Package

5 x 1.2ml syringes, disposable tips